



**Did you know
Metrodent manufacture
dental waxes ?**

No ? ...then take a look inside

Bite Blocks
Boxing In Wax
Gnathic Wax
Modelling Wax
Ribbon Wax
Soft Bite Wax
Strip Wax
Wax Up Sticks

Tried and Trusted



Metrodent Waxes

With over 80 years manufacturing experience Metrodent produce a comprehensive range of dental waxes tried and trusted by dentists and dental technicians both in the UK and around the world.

Using the finest raw materials Metrodent waxes are blended and produced to strict quality control standards giving you a first class British made product.



Made by Metrodent





Bite Blocks

For detailed occlusal records

Metrodent Bite Blocks can be easily adapted to shape and carved with a wax knife. Designed to represent the tooth position and are easy to use. Ideal for a wide range of temperature and climate conditions.

Colours	Pink Light Pink
Block Size	114 x 13 x 10 mm
Contents	48 Blocks (8 x 6)
Pack Size	585 g
Order Code - Pink	WAXBB48
Order Code - Light Pink	WAXBB48LP





Boxing In Wax

For boxing impressions prior to casting


Suitable for production of upper and lower jaws, for opposing models as well as individual split-cast bases. It surrounds the casting, can be easily adapted to its form and builds the edge for the model after the pouring of gypsum.

Colour	Red
Sheet Size	300 x 40 x 2.5 mm
Contents	18 Sheets
Pack Size	450 g
Order Code	BIW



Made by Metrodent

Modelling Wax

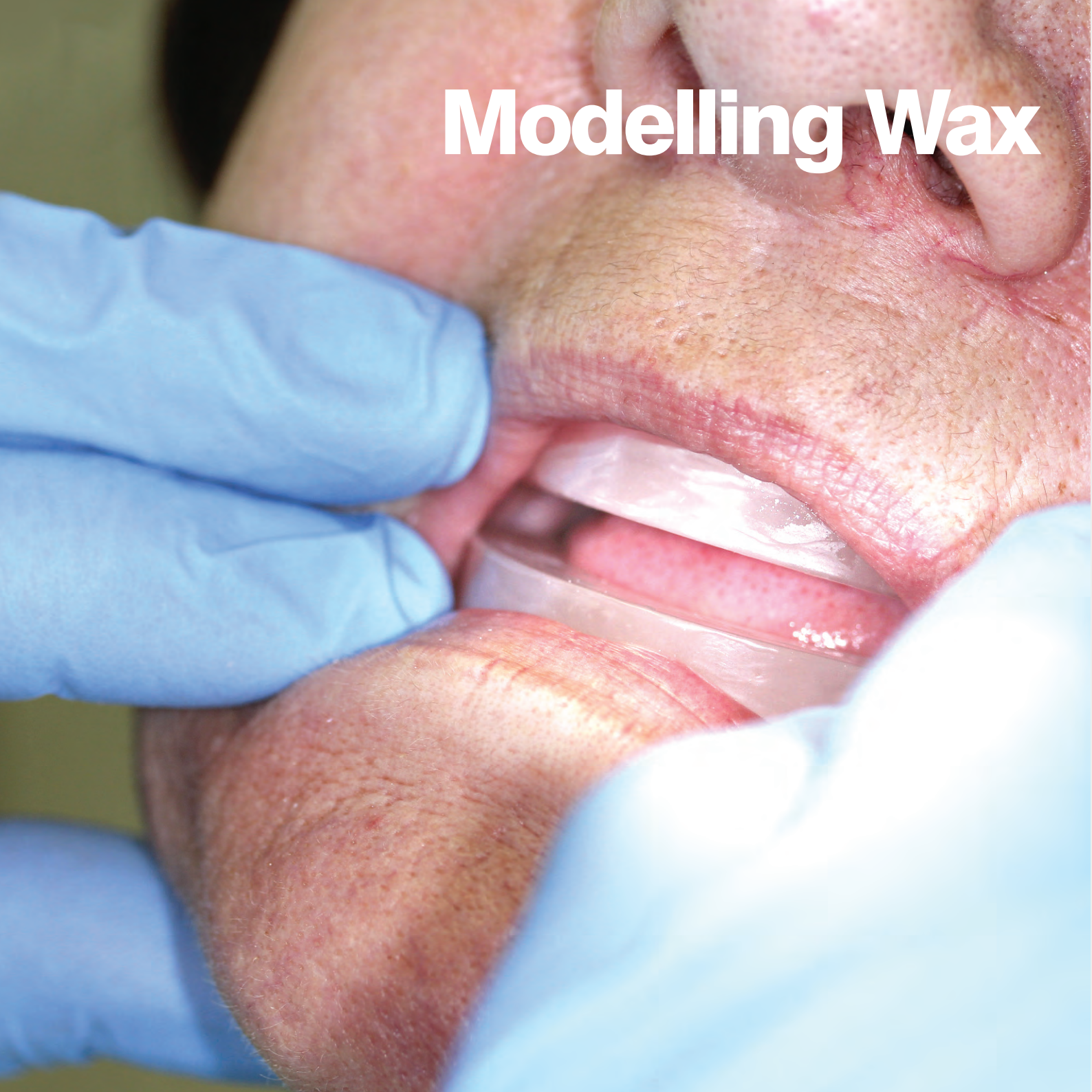

Gnathic

Differentiate your private work at the bite rim stage Available in Sheet Wax and Bite Blocks

Ivory wax sheets and preformed Ivory wax bite blocks for superior occlusal analysis.
More accurately represents the dentate area of the rim.
Improves patient communication and the intra oral contrast
makes the registration process easier.

Colour	Ivory
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Size	500 g
Order Code - 1.75 mm	GWAX500
Order Code - 1.50 mm	GWAX500EXPORT

Bite Blocks	114 x 13 x 10 mm
Contents	48 Blocks (8 x 6)
Pack Size	585 g
Order Code	GWAXBB48



Made by Metrodent

Modelling Wax



Light Pink

A very popular shade

The demand for Metrodent to produce a modelling wax in a light pink shade was overwhelming. Rapidly becoming a best seller our Light Pink modelling wax is based on our No.2 formula. A blend of paraffin and micro-crystalline waxes gives familiar ease of handling and consistent qualities aiding smooth work processing in the laboratory. This desirable colour option offers excellent soft tissue representation at try-in stage and is well received by both patients and clinicians.

Colour	Light Pink
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Sizes	500 g 2.5 Kg 20 Kg
Order Code - 1.75 mm x 500 g	LPWAX500
Order Code - 1.75 mm x 2.5 Kg	LPWAX2.5
Order Code - 1.75 mm x 20 Kg	LPWAX20
Order Code - 1.50 mm x 500 g	LPWAX500EXPORT
Order Code - 1.50 mm x 2.5 Kg	LPWAX2.5EXPORT
Order Code - 1.50 mm x 20 Kg	LPWAX20EXPORT



Modelling Wax



No.1

Highest melting point, ideal for carving

Our highest melting point wax originally designed for warmer climates. Utilising a higher melting point paraffin wax giving a longer molecular chain, it is more stable at higher ambient temperatures and therefore gives excellent carving properties. Ideally suited for high end prosthetics where the try-in is extensively contoured and detailed and where the try-in may be in the mouth for an extended period. Perfect for a hot laboratory in the summer months No.1 will be more brittle in lower temperatures.

Colour	Very Pale Pink
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Sizes	500 g 2.5 Kg 20 Kg
Order Code - 1.75 mm x 500 g	1WAX500
Order Code - 1.75 mm x 2.5 Kg	1WAX2.5
Order Code - 1.75 mm x 20 Kg	WAX20
Order Code - 1.50 mm x 500 g	1WAX500EXPORT
Order Code - 1.50 mm x 2.5 Kg	1WAX2.5EXPORT
Order Code - 1.50 mm x 20 Kg	1WAX20EXPORT



Made by Metrodent

Modelling Wax

A close-up photograph of a stack of several rectangular blocks of pink modelling wax. The blocks are stacked in a slightly offset manner, creating a sense of depth. The top block is in sharp focus, showing its smooth, slightly textured surface. The blocks below it are progressively more blurred, creating a strong sense of perspective. The lighting is soft and even, highlighting the uniform pink color of the wax.

No.2

Our most popular blend

This is our most popular modelling wax both in the UK and abroad. It is a blend of paraffin and micro-crystalline waxes and is designed to give good overall handling properties for everyday laboratory use. Pale Pink in colour No.2 is the longest manufactured blend in our range. Ideal for high workflow prosthetic laboratories where its consistency aids smooth production.

Colour

Pale 'Bluish' Pink

Sheet Sizes

183 x 87 x 1.75 mm

183 x 87 x 1.50 mm

183 x 87 x 3.50 mm

Pack Sizes

500 g

2.5 Kg

20 Kg

Order Code - 1.75 mm x 500 g

2WAX500

Order Code - 1.75 mm x 2.5 Kg

2WAX2.5

Order Code - 1.75 mm x 20 Kg

2WAX20

Order Code - 1.50 mm x 500 g

2WAX500EXPORT

Order Code - 1.50 mm x 2.5 Kg

2WAX2.5EXPORT

Order Code - 1.50 mm x 20 Kg

2WAX20EXPORT

Order Code - 3.50 mm x 500 g

DTW500

Order Code - 3.50 mm x 2.5 Kg

DTW2.5

Order Code - 3.50 mm x 20 Kg

DTW20



Made by Metrodent

Modelling Wax



No.4

Darker Pink variant on our most popular blend

This modelling wax is a colour option of our most popular wax.
Its handling properties and indications for use are the same as No.2.

Colour	Darker Pink
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Sizes	500 g 2.5 Kg 20 Kg
Order Code - 1.75 mm x 500 g	WAX500
Order Code - 1.75 mm x 2.5 Kg	4WAX2.5
Order Code - 1.75 mm x 20 Kg	4WAX20
Order Code - 1.50 mm x 500 g	4WAX500EXPORT
Order Code - 1.50 mm x 2.5 Kg	4WAX2.5EXPORT
Order Code - 1.50 mm x 20 Kg	4WAX20EXPORT



Made by Metrodent

Modelling Wax



Metromorphic

A pleasant aroma and easy to work

Metromorphic is a blend of paraffin wax, microcrystalline waxes and natural beeswax which with the omission of resins gives this wax a pleasant aroma. A useful everyday modelling wax.

Colour	Pale 'Bluish' Pink
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Sizes	500 g 2.5 Kg 20 Kg
Order Code - 1.75 mm x 500 g	7WAX500
Order Code - 1.75 mm x 2.5 Kg	7WAX2.5
Order Code - 1.75 mm x 20 Kg	7WAX20
Order Code - 1.50 mm x 500 g	7WAX500EXPORT
Order Code - 1.50 mm x 2.5 Kg	7WAX2.5EXPORT
Order Code - 1.50 mm x 20 Kg	7WAX20EXPORT



Made by Metrodent

Modelling Wax

The background of the slide is an abstract composition of soft, wavy lines in various shades of pink and red. In the upper-left quadrant, there is a rectangular area with a more pronounced, fine-grained texture, possibly representing a piece of wax or a specific material being modeled. The overall aesthetic is clean and modern, with a focus on organic, flowing shapes.

No.6

Easy to work with enhanced carving properties

Designed to be easy to work with No.6 is a blend of paraffin waxes, micro-crystalline waxes, natural beeswax and resins. The addition of resins enhances the carving properties of the wax making it suitable for contouring and characterisation of the denture base from try-in through to fit. Suitable for prosthetic laboratories with high workflow who may need to characterise their try-ins and require an all round reliable modelling wax.

Colour	Pale Pink
Sheet Sizes	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Sizes	500 g 2.5 Kg 20 Kg
Order Code - 1.75 mm x 500 g	6WAX500
Order Code - 1.75 mm x 2.5 Kg	6WAX2.5
Order Code - 1.75 mm x 20 Kg	6WAX20
Order Code - 1.50 mm x 500 g	6WAX500EXPORT
Order Code - 1.50 mm x 2.5 Kg	6WAX2.5EXPORT
Order Code - 1.50 mm x 20 Kg	6WAX20EXPORT



Made by Metrodent

Modelling Wax

College

For training and education purposes

Ideal for Dental Schools, Colleges and Universities this wax is made of secondary blends of our primary waxes and varies in colour and handling properties.

Colour

Pink

Sheet Sizes

183 x 87 x 1.75 mm

183 x 87 x 1.50 mm

Pack Sizes

500 g

20 Kg

Order Code - 1.75 mm/500 g

CWAX500

Order Code - 1.75 mm/20 Kg

CWAX20

Order Code - 1.50 mm/500 g

CWAX500EXPORT

Order Code - 1.50 mm/20 Kg

CWAX20EXPORT



Made by Metrodent



Ribbon Wax

For dorsal damming and other intra oral applications

Metrodent Ribbon Wax also known as Periphery Wax is used by dentists and technicians to extend the peripheral edge on an impression tray. High degree of plasticity and exact dorsal dam.

Colour	Red
Contents	126 Strips
Strip Size	210 x 8.0 x 1.50 mm
Pack Size	400g
Order Code	RWAX





Soft Bite Wax

Economical bite registration wax

Metrodent Soft Bite Wax is a proven bite registration material, formulated to provide dentists with a dense yet easily softened bite wax that is stable at room temperature. The thickness has been specially chosen to provide the most accurate registration.

Colour	Red
Contents	24 Strips (6 x 6 sheets)
Strip Size	230 x 12 x 7.5 mm
Pack Size	400g
Order Code	SRBW



Made by Metrodent



Strip Wax

An alternative bite registration wax

Metrodent Strip Wax is a proven bite registration material, formulated to provide dentists with a dense yet easily softened bite wax that is stable at room temperature.

The thickness has been specially chosen to provide the most accurate registration.

Colour	Red
Contents	35 Strips (7 x 5 sheets)
Strip Size	230 x 12 x 7.5 mm
Pack Size	200g
Order Code	SWAX



Made by Metrodent



Wax Up Sticks

For setting up and waxing try-in dentures

Metrodent Wax Up Sticks are specially formulated for the waxing up process of either the bite block or try in stages. It carves easily and boils out cleanly.

Colours

Pink

Contents

Light Pink

Strip Size

64 Sticks

Pack Size

ø 9 x 116 mm

Order Code - Pink

400g

Order Code - Light Pink

WUS

LPWUS



Made by Metrodent



Metrodent.com

