



## PractiPal®

Directa's products are developed, tested and evaluated by our selected team of highly-qualified and renowned Swedish and international dental professionals. The aim is to make life easier for all dental professionals in their daily clinical work.

**DIRECTA**





## **Suitable and practical: what are the ideal features of a tray system?**

Suitable and practical: what are the ideal characteristics of a tray system?

Infection control is a major challenge in healthcare practice. One area of high vulnerability is instrument reprocessing; therefore, a good tray system is essential. The success of cleaning, disinfection and sterilization is highly dependent on the design of the tray itself. A good tray system should facilitate the efficient handling of instruments throughout the complete reprocessing/cleaning cycle. An effective tray system with pre-prepared instrument sets used for each different type of treatment prevents the risk of confusing clean and used instruments, saves time and improves the clinic's efficiency.

Trays should be made of a light-non-metallic material. This avoids electrolytic corrosion, resulting from dissimilar metals being in contact with each other. It also reduces the heat energy loss and noise that can

occur from using metallic trays. The lightness of the tray reduces the time required for disinfection and sterilisation, as the heavier the load, the longer this process takes. Instruments should be held securely within the tray, but free to move for thorough cleaning. Trays that hold instruments too firmly, causing areas of contact, and/or have insufficient water penetration can leave microbiological residue on the instruments.

Very important and often overlooked are the logistics of instrument handling. A common mistake is overloading, which can make it difficult to keep track of instruments. PractiPal helps prevent this through its clever tray design, making it easier for you to be organized and maintain high standards of infection control. Regular reviews are recommended to ensure that cleaning efficiency, storage capacity and sterility are optimized. Overall, PractiPal is an ideal system to use in order to optimize your infection control process.







## **Enjoy the pleasure and efficiency of PractiPal® in your practice.**

Quiet, light-weight and colourful, PractiPal is a pleasure to work with! It is made of highly durable plastic that does not absorb heat in steam autoclaves and washer disinfectors, aiding the cleaning process and saving energy costs. PractiPal Tray System includes trays with bur stands, file stands, consumables and accessories that can easily be secured onto the trays throughout the entire cleaning and disinfection process. Both the tray clamp and bur stand/file holder are easy to open. Once prepared, a tray configuration can be cycled over and over again.

Our fibreglass-reinforced, durable plastic does not clang like metal can. Instruments are held securely in PractiPal, which can help avoid accidents while, at the same time, simplifying the inspection of instruments. Single-use materials and hazardous waste can be easily removed and thrown in the waste bin.

After treatment, PractiPal allows for the disposal of waste, but also holds instruments and burs in place for cleaning and sterilisation. Your pre-prepared tray contains only what you need for the procedure. Experience has shown that most dental practices have lots of superfluous instruments. Savings as high as 20% have been reported by dentists who have adopted this system and put all instruments in the practice into use.

Bur and accessory holders can easily be misplaced. Our simple solution is to secure the instruments, burs and accessories in pre-prepared trays for use over and over again. After treatment, the used but complete tray is replaced with a new configuration for the next patient.



## Organized work flow

Pre-configured tray systems minimize the quantity of instruments required whilst maximizing treatment efficiency, thanks to a simple cleaning procedure.

## Safety first

PractiPal holds all instruments securely, preventing accidents and cross-contamination.

## Improved hygiene control

The open design allows for maximum water and steam penetration, ensuring an optimal cleaning effect. Instruments are secure yet loose enough to avoid contact areas between the instruments and tray.

## Light and non-metallic

Avoids electrolytic corrosion. Lightness of the tray reduces heat energy loss and risk of overload in the autoclave.







**PractiPal®  
do not compromise on  
predictable cleaning.**

The belief that sterile is “clean and safe” is a common misconception. Residue can still remain on instruments even after cleaning, and this “sterile” residue can cause severe inflammatory reactions. Predictable cleaning is essential for 100% control. PractiPal trays, with a minimized frame, incorporate an ingenious means of securing instruments loosely. This ensures maximum water penetration for optimal washing, ultrasonic cleaning or simply rinsing.

Optimal cleaning, easy storage and safe handling makes PractiPal the product of choice. PractiPal includes bur, file and accessory holders that clip securely into the tray but are also easy to remove. The wide range of colours and accessories make it easy for you to design unique configurations to suit each surgery, dentist and individual operator. Once a cassette is configured, it can be reused time and time again.





PractiPal® Full Tray



## PractiPal® Full Tray

(280x188x37 mm)

Two half pieces assembled with joining clips. Organises and secures instruments. Holds 9 instruments. Bur/Endo stand holds 16 burs or 28 files. Available in nine colours. Colour of the tray: black. PractiPal Complete Set contains: Tray (2 pcs), clips (4 pcs), Instrument Clamp (1 pc), Bur/Endo stand for 16 burs or 28 files (1 pc), MultiBlock (1 pc), MultiDappen (5 pcs), MultiPoint (5 pcs), Waste Cup (5 pcs), Foam Large (1 pc) and Foam Small (5 pcs).

<b>Prod No</b>	<b>Product Description</b>
115199	PractiPal Complete Set, black/grey
115191	PractiPal Complete Set, black/red
115192	PractiPal Complete Set, black/yellow
115193	PractiPal Complete Set, black/blue
115194	PractiPal Complete Set, black/green
115195	PractiPal Complete Set, black/lilac
115196	PractiPal Complete Set, black/pink
115197	PractiPal Complete Set, black/mint green
115198	PractiPal Complete Set, black/orange





PractiPal® Half Tray



## PractiPal® Half Tray

(140x188x37 mm)

PractiPal Half Tray suits all set-ups and contains one half tray and one instrument clamp, holding up to 9 instruments. Available in nine colours. Colour of the tray: black

<b>Prod No</b>	<b>Product Description</b>
115100	PractiPal Tray, 2 clips
115106	PractiPal Tray and Instrument Clamp, black/lilac
115107	PractiPal Tray and Instrument Clamp, black/pink
115108	PractiPal Tray and Instrument Clamp, black/mint green
115109	PractiPal Tray and Instrument Clamp, black/orange
115110	PractiPal Tray and Instrument Clamp, black/grey
115111	PractiPal Tray and Instrument Clamp, black/red
115112	PractiPal Tray and Instrument Clamp, black/yellow
115113	PractiPal Tray and Instrument Clamp, black/blue
115114	PractiPal Tray and Instrument Clamp, black/green





PractiPal® Mini Tray



## PractiPal® Mini Tray

(80x188x37 mm)

PractiPal Mini Tray holds 5 instruments. This makes it ideal for check-ups, for dental hygienists and to stock spare instruments. Available in nine colours. Colour of the tray: black.

<b>Prod No</b>	<b>Product Description</b>
115020	PractiPal Mini Tray, Black
115050	PractiPal Mini Tray and Mini Clamp, black/grey
115051	PractiPal Mini Tray and Mini Clamp, black/red
115052	PractiPal Mini Tray and Mini Clamp, black/yellow
115053	PractiPal Mini Tray and Mini Clamp, black/blue
115054	PractiPal Mini Tray and Mini Clamp, black/green
115055	PractiPal Mini Tray and Mini Clamp, black/lilac
115056	PractiPal Mini Tray and Mini Clamp, black/pink
115057	PractiPal Mini Tray and Mini Clamp, black/mint green
115058	PractiPal Mini Tray and Mini Clamp, black/orange






## PractiPal® Instrument Clamp

Secures instruments for optimal cleaning without damage. Available in nine different colours and two sizes. One for Full and Half Tray and one for PractiPal Mini Tray.

Prod No	Product Description
115150	PractiPal Instrument Clamp, grey
115151	PractiPal Instrument Clamp, red
115152	PractiPal Instrument Clamp, yellow
115153	PractiPal Instrument Clamp, blue
115154	PractiPal Instrument Clamp, green
115155	PractiPal Instrument Clamp, lilac
115156	PractiPal Instrument Clamp, pink
115157	PractiPal Instrument Clamp, mint green
115158	PractiPal Instrument Clamp, orange

Prod No	Product Description
115040	PractiPal Mini Clamp, grey
115041	PractiPal Mini Clamp, red
115042	PractiPal Mini Clamp, yellow
115043	PractiPal Mini Clamp, blue
115044	PractiPal Mini Clamp, green
115045	PractiPal Mini Clamp, lilac
115046	PractiPal Mini Clamp, pink
115047	PractiPal Mini Clamp, mint green
115048	PractiPal Mini Clamp, orange

Select your  
colour coded  
treatment tray

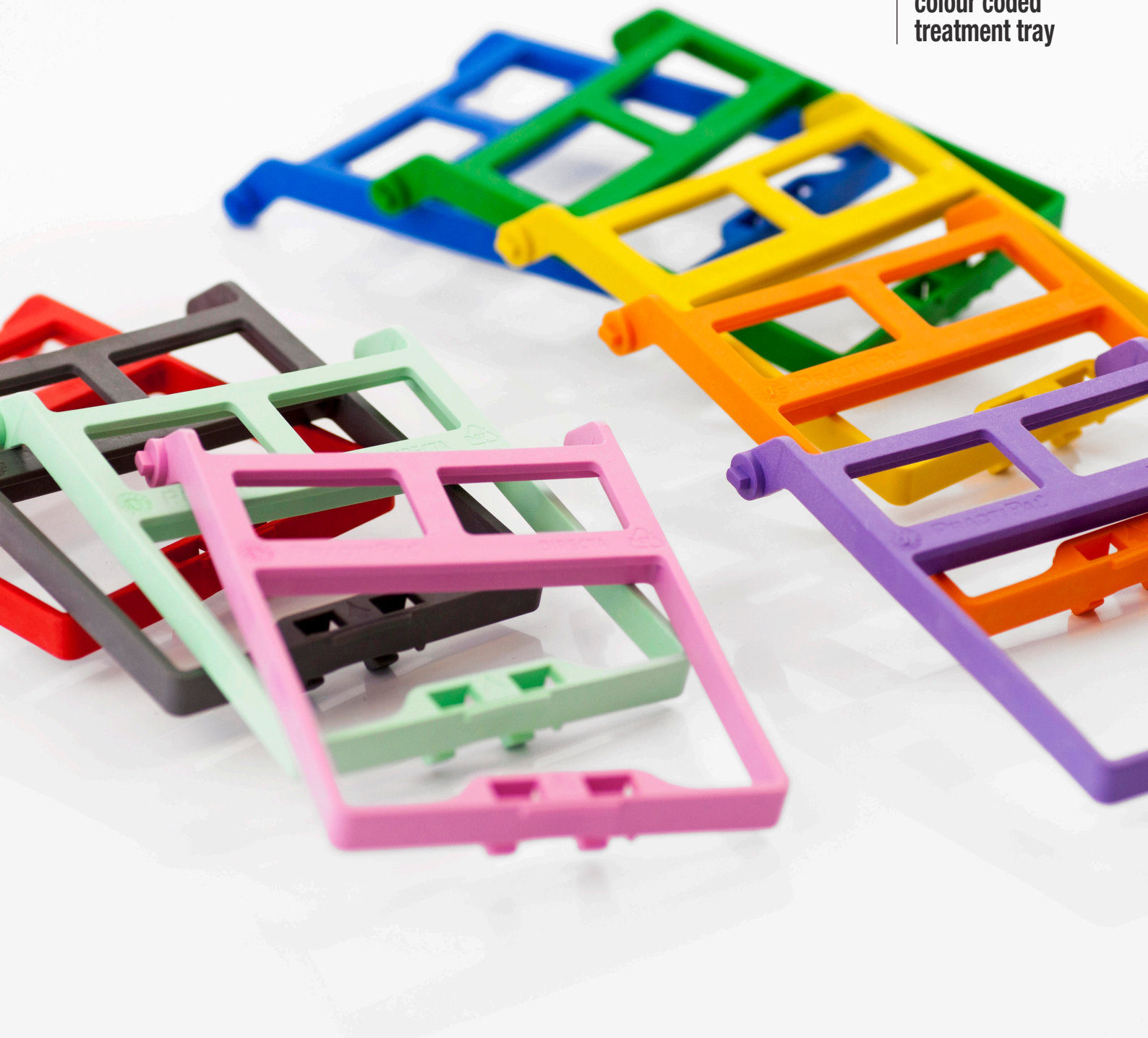





## PractiPal® Compact Bur Stand

Practical and safe bur holder for organization of burs.

Holds 2 long burs, 5 RA and 8 FG burs.

Easy handling - Minimizes sharps injuries.

Can be used together with the PractiPal system or as a stand alone.

- Prod No: 115121 PractiPal Compact Bur Stand, lilac, 1 pc
- Prod No: 115122 PractiPal Compact Bur Stand, pink, 1 pc
- Prod No: 115123 PractiPal Compact Bur Stand, mint green, 1 pc
- Prod No: 115124 PractiPal Compact Bur Stand, orange, 1 pc
- Prod No: 115185 PractiPal Compact Bur Stand, grey, 1 pc
- Prod No: 115186 PractiPal Compact Bur Stand, red, 1 pc
- Prod No: 115187 PractiPal Compact Bur Stand, yellow, 1 pc
- Prod No: 115188 PractiPal Compact Bur Stand blue, 1 pc
- Prod No: 115189 PractiPal Compact Bur Stand, green, 1 pc



## **PractiPal® Compact File Stand**

Compact file holder for the practical and safe organization of endodontic files.

The file stand holds 13 files.

Easy handling - Minimizes sharps injuries.

Can be used together with the PractiPal system or as a stand alone.

- Prod No: 115176 PractiPal Compact File Stand, lilac, 1 pc
- Prod No: 115177 PractiPal Compact File Stand, pink, 1 pc
- Prod No: 115178 PractiPal Compact File Stand, mint green, 1 pc
- Prod No: 115179 PractiPal Compact File Stand, orange, 1 pc
- Prod No: 115180 PractiPal Compact File Stand, grey, 1 pc
- Prod No: 115181 PractiPal Compact File Stand, red, 1 pc
- Prod No: 115182 PractiPal Compact File Stand, yellow, 1 pc
- Prod No: 115183 PractiPal Compact File Stand, blue, 1 pc
- Prod No: 115184 PractiPal Compact File Stand, green, 1 pc





## PractiPal® Bur/Endostand

Functioning as a bur stand or file holder, the same stand can hold either 28 files or 16 burs. Allows for easy inspection of files and burs, minimizing the risk of sharps injury.

Prod No	Product Description
115130	PractiPal Bur/Endostand, grey
115131	PractiPal Bur/Endostand, red
115132	PractiPal Bur/Endostand, yellow
115133	PractiPal Bur/Endostand, blue
115134	PractiPal Bur/Endostand, green
115135	PractiPal Bur/Endostand, lilac
115136	PractiPal Bur/Endostand, pink
115137	PractiPal Bur/Endostand, mint green
115138	PractiPal Bur/Endostand, orange



## PractiPal® MultiBlock

Holder that stores disposable items within its multiple inserts.  
Fits PractiPal Half and Full Tray.

PractiPal MultiBlock contain: PractiPal MultiBlock (1 pc), PractiPal MultiDappen (5 pcs), PractiPal MultiPoint (5 pcs), PractiPal Waste Cup (5 pcs) and PractiPal Foam small (5 pcs).

Prod No	Product Description
115140	PractiPal MultiBlock, grey
115141	PractiPal MultiBlock, red
115142	PractiPal MultiBlock, yellow
115143	PractiPal MultiBlock, blue
115144	PractiPal MultiBlock, green
115145	PractiPal MultiBlock, lilac
115146	PractiPal MultiBlock, pink
115147	PractiPal MultiBlock, mint green
115148	PractiPal MultiBlock, orange






## PractiPal® Half Tray Compact File Set

(140x188x37 mm)

Compact file holder stores 13 files. Capacity: 5 instruments. Available in nine colours.  
Colour of the tray: black PractiPal Compact File Set contains: PractiPal Half Tray,  
PractiPal Instrument Clamp for Half Tray, Compact File Stand and Waste Cups (4 pcs).

<b>Prod No</b>	<b>Product Description</b>
115200	PractiPal Compact File Set, black/grey
115201	PractiPal Compact File Set, black/red
115202	PractiPal Compact File Set, black/yellow
115203	PractiPal Compact File Set, black/blue
115204	PractiPal Compact File Set, black/green
115205	PractiPal Compact File Set, black/lilac
115206	PractiPal Compact File Set, black/pink
115207	PractiPal Compact File Set, black/mint green
115208	PractiPal Compact File Set, black/orange

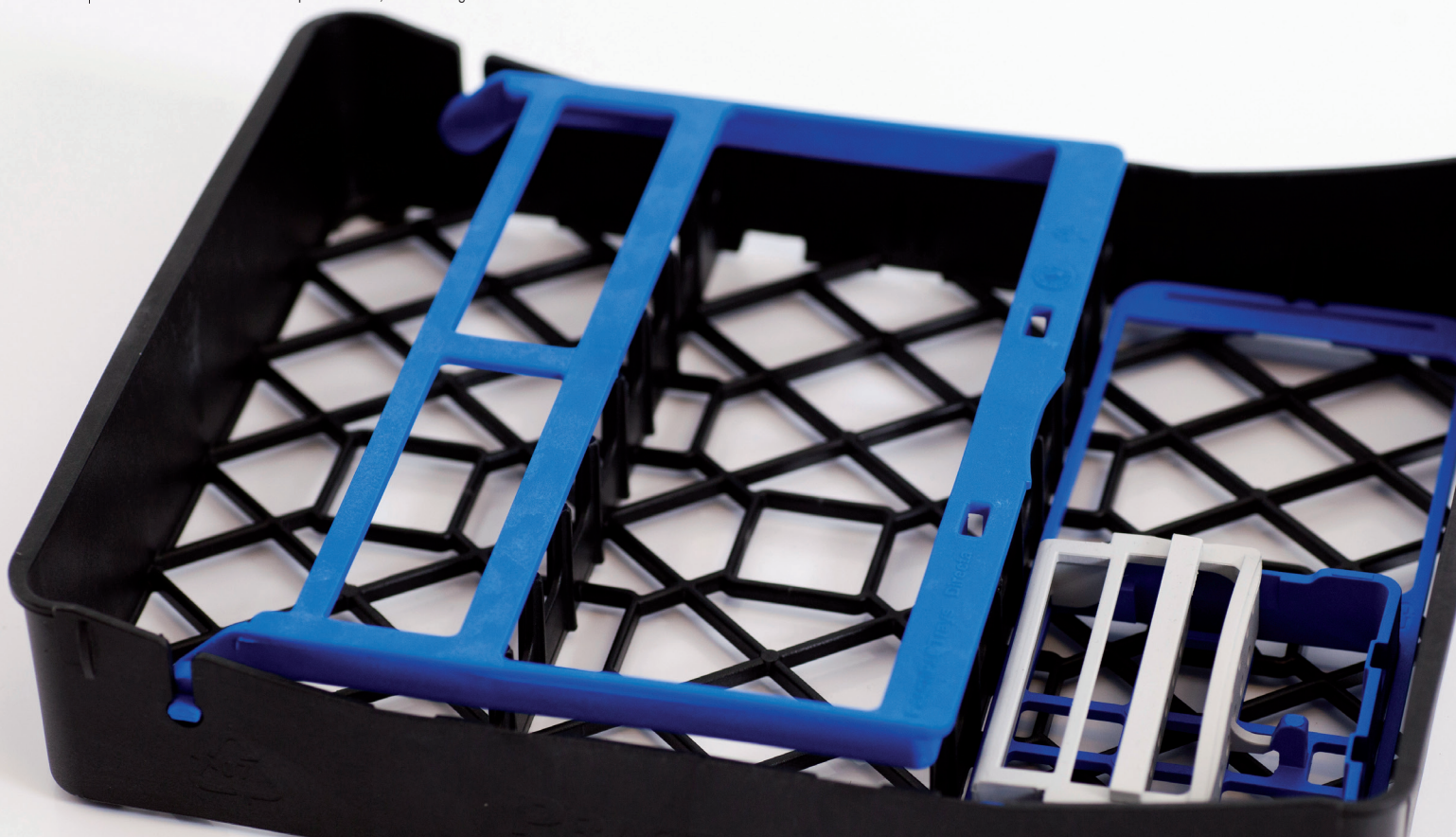



## PractiPal® Half Tray Compact Bur Set

(140x188x37 mm)

Compact Bur Set holds 2 extra long burs, 5HP/RA and 8FG. Capacity: 5 instruments. Available in nine colours. Colour of the tray: black PractiPal Compact Bur Set contains: PractiPal Half Tray, PractiPal Instrument Clamp for Half Tray, Compact Bur Stand and Waste Cups (4 pcs).

Prod No	Product Description
115160	PractiPal Compact Bur Set, black/red
115161	PractiPal Compact Bur Set, black/yellow
115162	PractiPal Compact Bur Set, black/green
115163	PractiPal Compact Bur Set, black/blue
115164	PractiPal Compact Bur Set, black/grey
115170	PractiPal Compact Bur Set, black/lilac
115171	PractiPal Compact Bur Set, black/pink
115172	PractiPal Compact Bur Set, black/mint green
115173	PractiPal Compact Bur Set, black/orange







## PractiPal® Waste Cup

Practical container for debris or liquids. Instruments may be scraped against the special jagged edges to remove debris. Can be used stand-alone or placed in the MultiBlock.

Prod No	Product Description
115002	PractiPal Waste Cup, 300 pcs



## PractiPal® Clips

Used to assemble two PractiPal Half Tray pieces together.

Prod No	Product Description
115120	PractiPal Clips, 2 pcs

## PractiPal® MultiDappen

Designed with indentations for different volumes of bonding agents, primers, etc. Can be used alone or placed on the MultiBlock.

Prod No	Product Description
115000	PractiPal MultiDappen, 300 pcs

## PractiPal® MultiPoint

Storage of Paper Points, Guttapercha Points and fluid. Can be used alone or placed in the MultiBlock.

Prod No	Product Description
115001	PractiPal MultiPoint, 300 pcs

## PractiPal® Foam

The unique slot design of the foam scraper allows for easy removal of debris from files and instruments. Large: 50x32x23 mm. Small: 30x15x13 mm.

Prod No	Product Description
115003	PractiPal Foam, small, 30x15x13 mm, 300 pcs
115004	PractiPal Foam, large, 50x32x23 mm, 50 pcs





PractiPal®

## Design By Dentists™

Directa's products are developed, tested and evaluated by our selected team of highly-qualified and renowned Swedish and international dental professionals. The aim is to make life easier for all dental professionals in their daily clinical work.

# DIRECTA

### DIRECTA AB

P.O. Box 723,  
194 27 Upplands Väsby, SWEDEN  
+46 8 506 505 75

### DIRECTA INC

64 Barnabas Road, Unit 3  
Newtown, CT 06470, United States  
203.491.2273

### DIRECTA GERMANY

Hammergutstr. 11  
92245 Kümmerbruck, GERMANY

info@directadental.com  
www.directadental.com  
www.directadentalgroup.com  
www.facebook.com/directadental  
www.instagram.com/directadental  
www.youtube.com/user/DirectaAB

DirectaDentalGroup

**DIRECTA** **topdental** **ORSING** **parkell** **physics**  
forceps

**RØNVIG** **TrollDental** **Sendoline** **CONTACTEZ**

2022-06\_REV01

